Silk'n®

Beautiful Technology

SONICSMILE (DELUXE)

SG-986




EN
DE
FR
NL
ES
PT
GR
TR
IT
NO
SV
DA

- EN QUICK INSTALLATION GUIDE
- DE INSTALLATIONSANLEITUNG
- FR GUIDE D'INSTALLATION
- NL INSTALLATIEHANDLEIDING
- ES GUÍA DE INSTALACIÓN
- PT GUIA DE INSTALAÇÃO RÁPIDA

- GR ΟΔΗΓΟΣ ΓΡΗΓΟΡΗΣ ΕΓΚΑΤΑΣΤΑΣΗΣ
- TR HIZLI KURULUM KILAVUZU
- IT GUIDA DI INSTALLAZIONE
- NO RASK INSTALLASJONSGUIDE
- SV SNABBINSTALLATIONSGUIDE
- DA HURTIG INSTALLATIONSGUIDE

Package contents/device parts

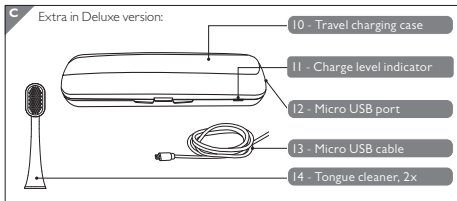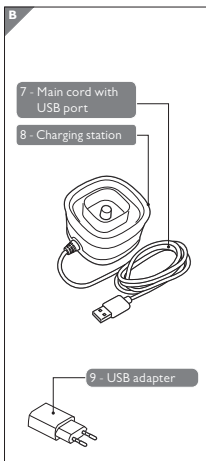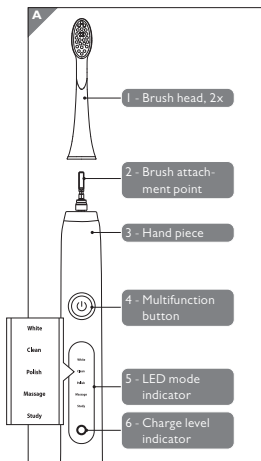*Note: The content of the box may vary based on the model purchased.

General information

These operating instructions accompany this SonicSmile (Deluxe) SG-986 (referred to below only as the “toothbrush”). It contains important information about start-up and handling. Before using the toothbrush, please read through the operating instructions carefully. This particularly applies to the safety instructions. You can also request these operating instructions in PDF format on our website located at www.silk'n.eu.

Safety

Only use the toothbrush as described in these operating instructions. Any other use is deemed improper and may result in damage to property or even personal injury. The toothbrush is not a children's toy. Silk'n accepts no liability for damage caused by improper or incorrect use.

First use

Clean all parts of the toothbrush before first use as described in the chapter “Cleaning”.

Operation

Charging the rechargeable battery

1. Insert the main cord with USB port ⑦ into the included USB adapter ⑨.
2. Place the toothbrush on the charging station ⑧. The toothbrush will now be inductively charged.
3. The charge level indicator ⑥ is illuminated in red during the charging process.
4. The toothbrush is fully charged after approx. 16 hours and the charge level indicator goes out. Once fully charged, the rechargeable battery will keep for up to 3 weeks.

Charging with the travel charging case (Deluxe version)

1. Place the hand piece ③ in the travel charging case ⑩. Insert the micro USB cable ⑬ into the micro USB port ⑫ of the case and the

other end into the USB adapter ⑨. Then, insert the USB adapter into a wall socket. The toothbrush will now be inductively charged. You can also charge via any CE certified device with a USB port. In this case, insert the micro USB cable into the micro USB port of the case and the other end into the USB port of the device.

2. The charge level indicators ⑥ ⑪ on the toothbrush and case are illuminated in red during the charging process.
3. The toothbrush is fully charged after approx. 16 hours and the charge level indicators go out. Once fully charged, the rechargeable battery will keep for up to 3 weeks.

Using the toothbrush

- Do not apply too much pressure to the toothbrush and your teeth when brushing them.
- Your gums may bleed the first few times you use the toothbrush. If your gums bleed for more than 14 days, please consult a dentist.
- Brush your teeth using small circular movements.
- Clean up to the spaces between your teeth.

Brushing your teeth

1. Attach a brush head ① to the hand piece ③ so that the brush attachment point ② of the hand piece extends into the brush head.
2. Select the desired mode by pushing the multifunction button within 3 seconds (see chapter "Selecting the mode"). Wait for more than 3 seconds to switch the toothbrush off.
3. Brush the upper and lower parts of your teeth as well as the exterior and interior of your teeth for 30 seconds each without applying pressure.

After every 30 second interval, the toothbrush will stop for 0.5 seconds and automatically switch off after 2 minutes.

4. After each use, clean the toothbrush as described in the chapter "Cleaning".

Brushing your tongue with the tongue cleaner (Deluxe version)

1. Wet the tongue cleaner ④ and brush your tongue with the micro bristle side in a back and forth motion for 30 seconds. Clean as far back as comfortable.
2. Rinse your mouth and use the scraping side on your tongue for 30 seconds.

Selecting the mode

The LED mode indicator ⑤ displays the current mode.

Function	Operation	Description
White	Push 1×	For white teeth and less staining.
Clean	Push 2×	For removing plaque, perfect for sensitive teeth and gums.
Polish	Push 3×	For polishing teeth.
Massage	Push 4×	For gently massaging your gums and improving the health of your gums.
Study	Push 5×	Setting for getting used to the toothbrush. The vibration intensity gradually increases from weak to strong during the first 12 applications.

Cleaning

1. After every use, hold the brush head ① under running water and switch the toothbrush on.
2. Remove all leftover toothpaste from the brush head and the hand piece ③ and switch the toothbrush back off.
3. Pull the brush head off of the hand piece.
4. Also clean the brush head by holding it under running water.
5. Dry the brush head and the hand piece. Clean the charging base and/or charging case with a soft, damp cloth and wipe it dry.

6. Change the brush head no later than every 3 months.

As soon as the colour of the brush head fades, change it. You can purchase new brush heads on our website at www.silkn.eu.

Storage

- Always store the toothbrush in a dry area.
- Use the original packaging, a similar box or the travel charging case for storage.

Troubleshooting

Problem	Possible cause	Solution
The toothbrush isn't charging.	The charging station or charging case is defective.	Try other sockets or USB connector jacks and charge the toothbrush for at least 16 hours.
	The USB connector jack is defective.	Try other sockets or USB connector jacks and charge the toothbrush for at least 16 hours.
	The charging station or charging case is defective.	Replace the charging station or charging case within the warranty period.
The toothbrush operates with only a small amount of vibration.	The rechargeable battery is depleted.	Charge the toothbrush for at least 16 hours.

Problem	Possible cause	Solution
The toothbrush does not run for 2 minutes.	The rechargeable battery is depleted.	Charge the toothbrush for at least 16 hours and make sure that the toothbrush is properly placed on the charging station or in the charging case.

Technical data

Model number: SG-986

Charging station

Input: DC 5V
 Power consumption: 1.5W
 Protection class: II
 IP code: IPX7
 Dimensions (L × W × H): 6.2 × 6.2 × 4.3 cm

Hand piece

Rechargeable battery: 3.7V, 700 mAh
 Operating time: 42 times/full charge
 Timer: 0.5 second pause after approx.
 30 seconds, deactivation after 2 minutes
 Dimensions (L × W × H): 2.6 × 2.7 × 25.4 cm

Adapter

Model number: HNFG050060UE
 Input: 100–240V~, 50/60 Hz, 0.2A
 Output: 5.0V \equiv 600mA
 Protection class: II

Travel charging case (Deluxe version)

Model number: SG-209
 Working voltage: DC 5V \equiv 300mA

Power rating:	1.5W
Dimensions (L x W x H):	22.6 x 6.4 x 3.6 cm

Warranty

This product is covered by a 2-year warranty according to European regulations and law. The extent of the warranty on this product is limited to technical defects caused by faulty production processes. In the event you want to claim warranty please be sure to contact our customer service for instructions. They may be able to solve your problem without needing to return the product to the store or our service center. Our customer service is always pleased to help you!

Customer service

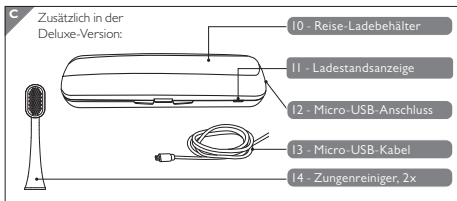
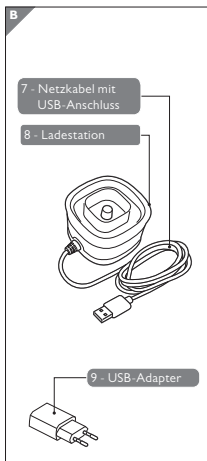
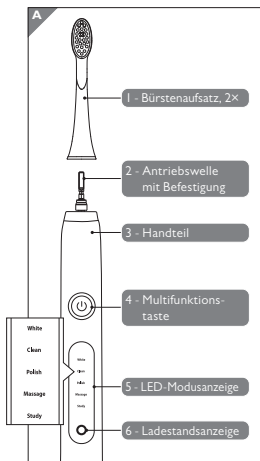
For more information about the products of Silk'n visit your regional Silk'n website: www.silkn.eu. If the device is damaged, defective, in need of repair or you need our assistance, please contact your nearest Silk'n service center. The full manual is also available as a PDF download from: www.silkn.eu.

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Invention Works B.V.
InnoEssentials Deutschland GmbH