
BODYSCALE

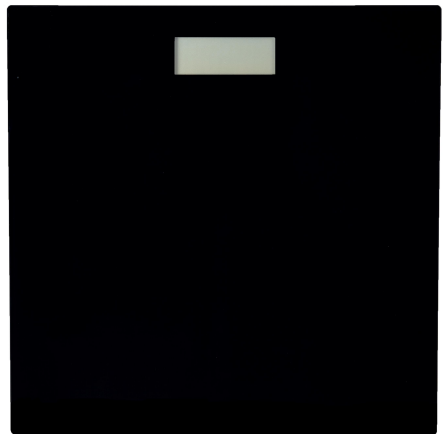Item No. 19192

Instruction manual

Please read all instructions carefully and retain for future use!

(EN) Content

- 1 x Scale
- 1 x Manual

Specifications

Dimensions	28 x 28 x 2.5 cm
Material	Glass + PS
Unit	kg/lb
Power input (DC)	3 V
Batteries	2 x AAA
Capacity	5kg - 180kg
Graduation	0.1 kg/0.2 lb
Low battery display	Lo
Overweight display	oL

Instructions

- Turn on auto
- Auto off in 10 seconds
- Remove the insulation film from the back of the scale near the battery compartment prior to operation.

- Place the scale on a plane and hard surface prior to weighting
- Step on the scale, balance your weight evenly and wait until the display reads your weight. Your weight will flash on the readout for 10 seconds.
- Press the button near the battery box and you can choose the different unit: kg/lb
- Replace the battery immediately when you see the Lo battery indicator.

ATTENTION:

- Do not leave any object on the scale when not in use.
- Keep the scale surface and your feet drying when weighting.
- Make sure you stand in the center of the scale to avoid tipping.
- Avoid dropping or hitting the scale.
- Wipe scale surface with a clean dry cloth if need clean.
- Do not wash the scale with water, the scales mechanism can be damaged if immersed in water.

WEEE:

This product bears the selective sorting symbol for Waste Electrical and Electronic Equipment (WEEE). This means that this product must be handled pursuant to European Directive 2012/19/EU in order to be recycled or dismantled to minimize its impact on the environment. The user has the choice to give this product to a competent recycling organization or to a retailer when buying a new electrical or electronical piece of equipment.

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