


WELCOME

Congratulations on your new Hypervolt Plus

You are one step closer to moving better,
feeling better, and being a better version of you.
The Hypervolt Plus was created to keep you
doing what you love. Whether you're an elite
athlete, weekend warrior, or simply trying to get
your body to feel its best, the Hypervolt Plus will
help you get there... and then some.

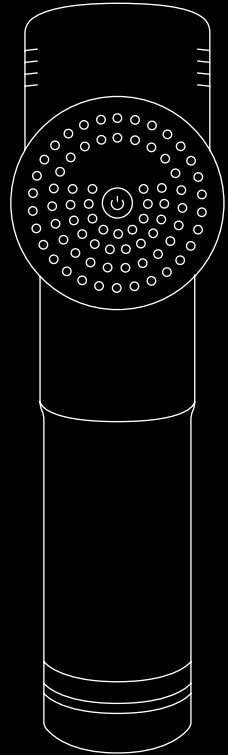

NEXT STEPS

Let's get started

1 Register your device at hyperice.com/register-product to activate your warranty and ensure easy returns, repairs, or refunds.

2 Download the **Hyperice App** to get the most out of the Hypervolt Plus, or any Hyperice device.



3 Charge your device by plugging in your new Hypervolt Plus with the provided Hyperice charger. Fully charge the battery for up to six hours before first use.

Have questions? Call 1-855-734-7224, chat live at hyperice.com, or email customersupport@hyperice.com.


QUICK START GUIDE

Moving better, made simple

1

Power up

Switch your Hypervolt Plus to the ON position at the base of the handle. Once your device is powered ON a blue light will illuminate, plus (up to) four green lights for battery level.

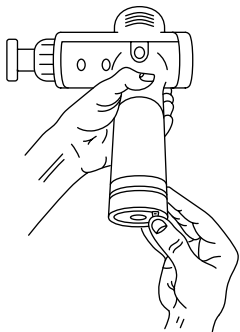

2

Activate pressure sensor

Activate the pressure sensor by pressing the speed setting button once. Wait five seconds as the white lights sequence and turn off.

Your Hypervolt Plus is now ready for use.




3

Adjust speeds


Choose your desired speed (1-3) by pressing the speed setting button once per speed.

The blue lights above the speed setting button will indicate the current speed.

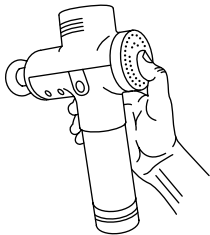**4**


Changing attachments

Choose desired head attachment and insert directly into the shaft / opening while pushing firmly. **Fork attachment?** Align the "groove" on the fork with the recess on the shaft / opening and push firmly.

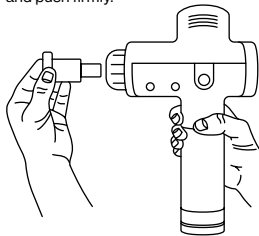**5**


Recharging

Only use the provided Hyperice wall charger to uphold warranty. Full charge is indicated when the LED "band" on the handle stops flashing and becomes solid green.

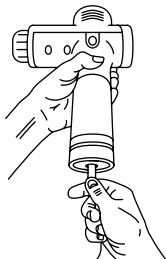**6**


Caring for your Hypervolt

Use a damp, clean, cloth and gently wipe your Hypervolt Plus. Make sure the power is OFF and the battery charger is not attached. Store in a clean, cool place when not in use.

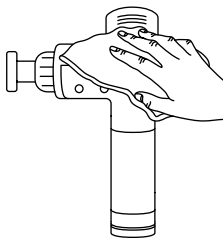


HEAD ATTACHMENTS

Choosing the best head attachment

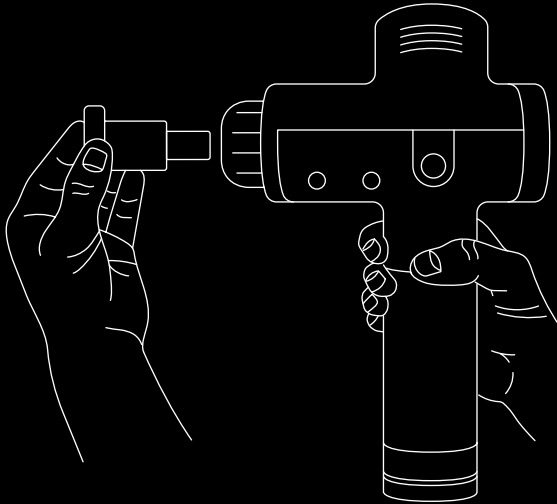

HEAD ATTACHMENTS

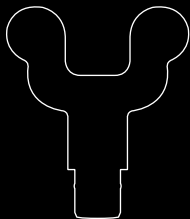


Fork

Longer muscles such as calves, forearms, biceps, and ankle / Achilles area.




Bullet

Trigger points found in glutes, hips, traps, shoulders, and other complex muscle groups.

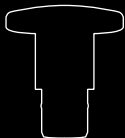


Flat

Larger surface areas such as quads, hamstrings, back, and chest.

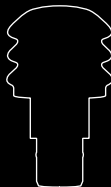

Cushion

Sensitive areas such as wrists, ankles, knees, traps, and around the neck.

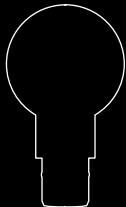

Ball

Sore muscles and a great starter attachment for larger muscle groups.

Want to learn more? Visit hyperice.com/how-to for video tutorials on your Hypervolt Plus or any Hyperice device.

HYPERICE APP

Your digitally connected wellness guide

Connect, monitor progress, learn,
and operate Hyperice products
from your IOS devices.



Download on the
App Store

HYPERICE APP

A few of the features



Apple Health™ and Strava™ connected

Connect to like-minded apps for a more personalized digital experience.



Data tracking / analytics

Tap into activity-based analytics to monitor your progress.



Personalized programs

Guided warmup, recovery, and maintenance sessions.



Operate your devices

Operate your Hypervolt or any Hyperice Bluetooth® connected device.



IOS is a registered trademark of Cisco and is used under license. The Apple logo, Apple Health, and the App Store are registered trademarks of Apple, Inc. and used under license. Strava is a registered trademark of Strava, Inc. and is used under license. The Bluetooth® word mark and logos are a registered trademark owned by Bluetooth SIG, Inc. and is used under license.

HYPERCARE

A dedicated team just for you

Every Hyperice product is backed by a 1-year warranty and award-winning support from our Hypercare™ team – a group of specialists dedicated to your overall well-being and expert guidance on Hyperice products.



HYPERCARE

Here's a few of the things your team can do for you




Product knowledge

Answer questions about your device and best practices.



Place new orders

Help with placing orders online and/or authorized resellers.

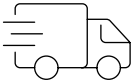

Track new orders

Track new orders and work with carriers to find solutions.




Billing solutions

Research and provide solutions for any billing issues.



Trouble shoot

Help trouble shoot and solve issues with your device.

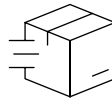

Initiate returns

Provide instructions for returns, refunds, or repairs.

Get in touch.

Hypercare® is registered trademark of Hyperice, Inc. Copyright 2020 © Hyperice, Inc. All rights reserved.

Call 1-855-734-7224

Chat [hyperice.com](https://www.hyperice.com)

Email customersupport@hyperice.com

Hyperice, Inc

525 Technology Drive, Suite 100


Irvine, California 92618

OTHER PRODUCTS

A family of products designed to upgrade the human body


Hypervolt

Percussion Massage Device

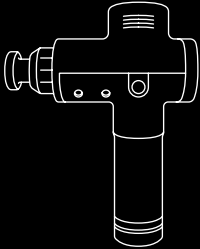


Vyper 2.0

Vibrating Fitness Roller

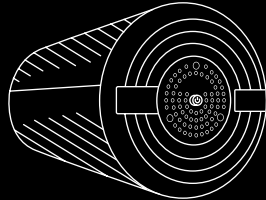


Hypersphere

Vibrating Massage Ball

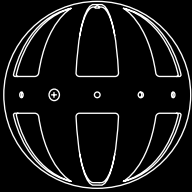


Hypersphere Mini

Compact Vibrating Massage Ball

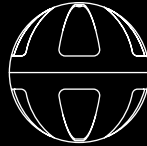


Venom

Wearable Heat + Vibration Devices

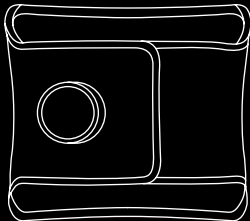

ICT

Ice Compression Devices

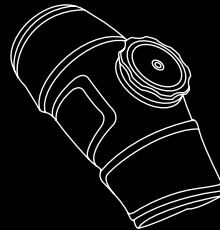

Learn more / shop @ hyperice.com