


# ASSAMBLY INSTRUCTIONS

## Sit-stand assistance serie 235

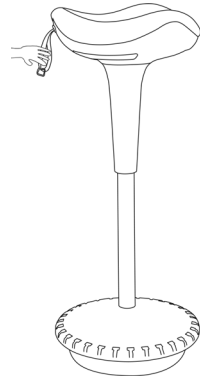

**WARNING:** FAILURE TO COMPLY WITH OR OBSERVE ALL ASSEMBLY, SAFETY AND OPERATION INSTRUCTIONS AND WARNINGS REGARDING THE USE OF THIS PRODUCT MAY RESULT IN SERIOUS BODILY INJURY.

### Important information

Read the entire instruction manual before you can start installation and assembly.

This product is not intended for use by young children without supervision.

This product is intended for indoor use only. Using this product outdoors could lead to product failure and personal injury.


**Important:** Ensure that you have received all parts according to the component checklist prior to installation. If any parts are missing or faulty, telephone your local distributor for a replacement.

Do not stand on this chair, do not use this chair as a step ladder.

Only use this product for placing one person at a time.


Do not use this chair unless all bolts, screws and knobs are properly tightened. Failure to heed these warnings could result in serious injury.

### COMPONENTS

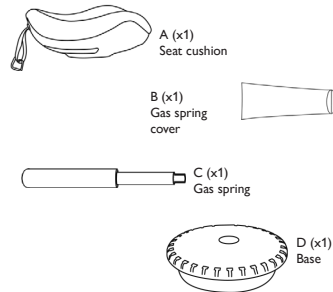


### 1- Assemble seat cushion

- Turn the seat cushion (A) upside down.
- Press the gas lift cover (B) into the holes on the bottom of the seat cushion.
- Turn counterclockwise to attach the gas lift hood (B) with seat cushion (A).

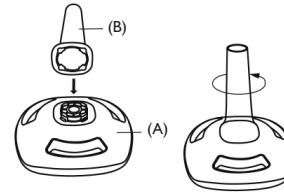


### 2- Assemble gas spring

- Remove the protective sleeve before mounting the seat cushion.

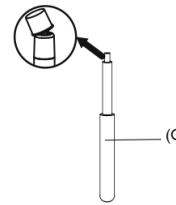


### 3- Assemble base

- Insert the gas lift (C) into the center hole of the base (D).

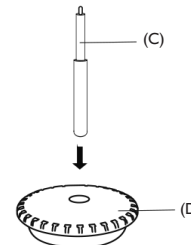

### 4- Final assembly operation

- Place the seat cushion (A) over the gas lift (C) and push it down firmly until the connection is secure.
- Your sit-stand assistance is now complete and ready to use

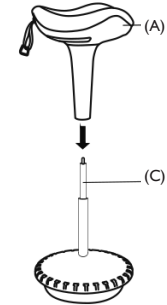

### 5- Adjust seat height

- While seated, lift your body weight off the chair. Slowly pull the strap up so that the crank reaches the desired height.
- To lower the seat height, slowly pull on the strap while seated.

